




## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification <sup>6</sup> : G01N 33/86, 33/96	A1	(11) International Publication Number: WO 95/30154 (43) International Publication Date: 9 November 1995 (09.11.95)
(21) International Application Number: PCT/US95/05195 (22) International Filing Date: 28 April 1995 (28.04.95) (30) Priority Data: 08/235,016            28 April 1994 (28.04.94)            US (71) Applicant: DADE INTERNATIONAL INC. [US/US]; 1717 Deerfield Road, P.O. Box 778, Deerfield, IL 60015-0778 (US). (72) Inventors: WOODHAMS, Barry, J.; 20, route Henri-Dunant, CH-1700 Fribourg (CH). BURGESS-WILSON, Michael, E.; 39, route du Coteau, CH-1752 Villars-sur-Glane (CH). (74) Agents: PEARSON, Louise, S. et al.; 1717 Deerfield Road, P.O. Box 778, Deerfield, IL 60015-0778 (US).	(81) Designated States: AU, CA, JP, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE). Published With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.	

(54) Title: CALIBRATOR FOR PROTHROMBIN TIME (PT) ASSAYS

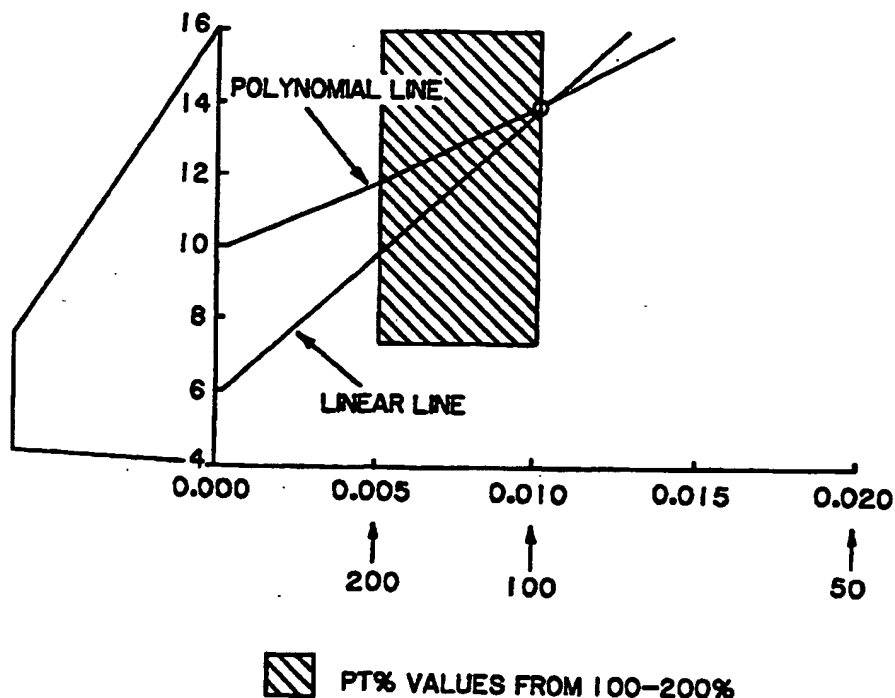

## (57) Abstract

This invention pertains to a PT Assay Calibrator and a method of preparing a PT Assay Calibrator including a coagulation factor such as recombinant FVII or recombinant FVIIa that will allow preparation of PT calibration curves with values about 100 % and which will give results analogous to those obtained using fresh normal plasma.

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LK	Sri Lanka	SN	Senegal
CN	China	LU	Luxembourg	TD	Chad
CS	Czechoslovakia	LV	Latvia	TG	Togo
CZ	Czech Republic	MC	Monaco	TJ	Tajikistan
DE	Germany	MD	Republic of Moldova	TT	Trinidad and Tobago
DK	Denmark	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	US	United States of America
FI	Finland	MN	Mongolia	UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

1

## CALIBRATOR FOR PROTHROMBIN TIME (PT) ASSAYS

Background of the InventionField of the Invention

This invention relates to a method of preparing a  
5 commercial plasma preparation that will allow  
preparation of PT calibration curves with values about  
100% and which will give results analogous to those  
obtained using fresh normal pooled plasma.

Description of the Related Art

10 The Prothrombin Time (PT) is used as a screening  
test for blood coagulation factor deficiencies and for  
monitoring oral anti-coagulant therapy using, e.g.,  
coumadin. Thromboplastin reagents activate the  
"extrinsic" pathway of coagulation and are the basis  
15 for the PT test. Thromboplastin contains lipidated  
tissue factor (TF), which is the activator of the  
extrinsic pathway. This activation centers on Factor  
VII (FVII) and activated Factor VII (Factor VIIa), the  
TF-FVII Complex activates Factor X, which with Factor  
20 V activates Factor II to produce thrombin, which  
creates the fibrin clot.

There are several ways of expressing the results  
of the PT test. One system, the INR system, is  
recommended by the World Health Organization.  
25 However, many countries have not adopted this system  
for expressing PT results. Moreover, the INR system

has only been validated for patients on oral anticoagulant control, but should not be used in expressing results from patients with other disease states, such as liver disease. Another system, commonly used in the United States, expresses the time in seconds for the blood to begin to coagulate. Still another system expresses the results in terms of a percentage PT ("% PT") which is read from a standard calibration (or dilution) curve prepared by diluting fresh normal pool plasma ("FNP") in 0.9% saline. (Other diluents work, but by convention, only saline is used.) The curve allows for the conversion of results from time in seconds to percent of normal activity (% PT). Unfortunately, in order for this system to be used, most laboratories have to prepare their own pool plasma and keep it frozen, usually in liquid nitrogen or frozen at -80°C. Moreover, due to the inherent variation found in different plasma pools, there is no standardization between the plasma pools of different laboratories. Moreover, it has been shown that if a pool of plasma is prepared, the mean % PT value obtained from the pool is different than the mean % PT value obtained from the individual samples that were used to make the pool. It has also been shown that the collection of bulk collections of blood, as would be required to commercially prepare a

lyophilized standard, causes a reduction in the measured % PT when compared with blood collected by venipuncture. See Important Differences Encountered in the Normal Plasma Pools used for the Control of Oral Anticoagulation. M. Burgess-Wilson, R. Burri and B. Woodhams, Thromb. Haemost. 69 Abstract 2081 (1993). Moreover, FNP cannot be sold until lyophilized. Lyophilization results in a plasma which, when reconstituted, has a % PT value lower than that found in a normal hospital pool of plasma. This reconstituted FNP is then used to prepare the standard curve. The dilutions usually used for the standard calibration curve are undiluted, 1:1, 1:2 and 1:4. Where reconstituted FNP is used, the undiluted sample is assigned a value of 100% PT. A PT assay is performed and the results (in seconds) are plotted on hyperbolic or reciprocal graph paper against the dilution (in %). See Fig. 1. Patient samples are tested undiluted and then read from this standard curve. However, using reconstituted FNP as a calibrator means that values for normal samples are above the top calibration point of the standard curve made using the reconstituted FNP. (By definition, 50% of all normal values would be above the top point of the standard curve.)

The % PT curve is not a straight line. Although a polynomial plot gives the most realistic curve through the data, many laboratories and users do not have the computer software required for such a procedure. Therefore, a linear curve through the points is commonly used. To make the results more accurate around the 100% region of the curve, the line is forced through the 100% point. One type of assay machine, the Medical Laboratory Automation ("MLA") Electra automated coagulometers, does not calculate % PT outside of certain ranges (above about 125% PT and below about 12% PT). The recommended method of calculating % PT varies between the instrument manufacturers. There is no universally used standard procedure. Some instrument manufacturers, such as MLA, recommend forcing a straight line through 100%. Others recommend polynomial or non-forced straight lines. This introduces variability into the procedure, especially if the calibration plasma has a value of % PT much lower than 100%. See Fig. 2. In the examples that follow, the method of calculating the %PT was to use a forced linear curve through the 100% point using the SigmaPlot transformation.

Summary of the Present Invention

This invention relates to a method for preparing a commercial plasma preparation that will allow calibration curves to be prepared that will have % PT values of about 100% and will give results analogous to those obtained with FNP. In summary, the invention involves the addition of recombinant human FVII or FVIIa (or any other source of FVII, provided it is of high enough purity and behaves in a similar fashion to human FVII) to normal human citrated plasma to give the required PT%. For an article discussing the purification of recombinant human Factor VII, please see Kemball - Cook, McVey, Garner, Martin, O'Brien and Tuddenham, Stable High Level Expression of Recombinant Human Factor VII In Mammalian Cell Culture, Thromb. Haemostatis 69 (6) 1993, Abstract 253. The resulting plasma is lyophilized and calibrated. It is expected that the addition of other recombinant factors such as rFVIII, rFV or rFXI could be made to a plasma that would also act as a calibrator for other coagulation assays, e.g. FVIII, FV, FXI, derived fibrinogen, FIX, FII, FX, Protein-C, Protein-S, and APTT (clotting and chromogenic) assays. For example, rFV is obtained from available sources and can be added to the plasma such that a level of about 100% rFV is achieved.

The calibrator of the present invention can be used with thromboplastin reagents such as Thromboplastin IS (Baxter's dried rabbit brain with calcium PT reagent), Thromboplastin C, and 5 Thromboplastin C+. Particularly, the calibrator of the present invention is designed for use with a recombinant tissue factor PT reagent such as Baxter Diagnostics Inc.'s Dade Innovin™ dried recombinant human tissue factor with calcium and Ortho Diagnostics 10 Systems Ortho® RecombiPlastin™ recombinant tissue factor relipidated with highly purified phospholipids reagent which are used as reagents in the PT determinations and PT-based assays. Recombinant tissue factor reagents, and in particular, Innovin™ 15 reagent, was found to have increased sensitivity, when compared to other reagents used in PT determination and PT-based assays, to various factor deficiencies and oral anticoagulant-treated patient samples. The increased sensitivity of such reagents is such that 20 they differentiate much more between FNP collected by syringe or by blood bag than traditional thromboplastins (prepared from animal or human tissue extracts). A calibration plasma should be collected in a fashion similar to clinical samples, i.e., 25 syringe drawn. However, until the calibrator of the present invention, commercial preparations of a



7

calibration plasma with a PT of 100% were difficult, if not impossible, to prepare. Lyophilized normal plasma has a %PT of 85% or less when measured with Innovin<sup>™</sup> reagent. The use of a plasma sample with

5 such a low % PT value makes calculations of the % PT value of normal samples difficult and introduces a large amount of variation according to the method used to calculate the % PT, as explained further below. As shown in Fig. 2, the boxed area shows the two curves

10 which can be drawn (polynomial and extrapolated). The enlarged boxed area shown in Fig. 3 demonstrates that the two curves will give very different results as they diverge. The divergence increases above the top calibration point. If the top calibration point is at

15 85%, then the calibration of normal results (130-70% PT) will be more strongly influenced by the choice of curve as the 85-100% PT part of the curve will have to be extrapolated. The use of the calibrator of the present invention will keep the % PT close to 100% and

20 avoid using the diverging areas of the curves. The resulting calibrated plasma preparation can be used on the MLA Electra, KC, and ACL range of instruments.

#### Detailed Description of the Drawings

Fig. 1 depicts a calibration curve of clotting

25 time (in seconds) vs. 1/PT% for dilutions of FNP in saline.

Fig. 2 depicts the problem posed by calculating % PT using PT dilution curves when the clotting time of the test plasma is shorter than that of the calibration plasma.

5 Fig. 3 depicts an enlarged portion of Fig. 2.

Fig. 4 depicts the effect of the addition of rFVIIa in different concentrations on PT Clotting Time in seconds from the data in Table 1a.

Fig. 5 depicts the PT calibration curves of FNP  
10 alone and FNP with the addition of rFVIIa ( $1/10^3$  dilution), from the data in Table 1b.

#### Detailed Description of the Invention

Recombinant FVIIa did raise the % PT of the plasma pool. Recombinant FVII also raised the % PT.  
15 The amount of recombinant material needed to be added to a large pool of plasma to produce a % PT of about 100% was determined.

The FVII levels achieved (as measured using the one stage clotting assay) did not usually parallel the  
20 rise in PT%. Two lots of rFVII showed quite different relationships between PT% and FVII level rise. The difference was thought to be due to "contamination" of the rFVII with the more active rFVIIa. As described later herein, the rFVIIa material did not have this  
25 problem. Without limiting the scope of the invention, it is believed that rFVIIa is preferable as a

calibrator because the "contamination" factor is not present.

Either rFVII or rFVIIa was added to a pool of HEPES buffered citrated plasma. While HEPES buffer was chosen for these examples because it lyophilizes well, it is believed that most buffers which work in the physiological pH range could be used, except for phosphate type buffers. Examples of buffers which should work include Good's Buffers: PIPES, ACES, BES, MOPS, TES, and TRICINE. The resulting plasma plus recombinant material was tested for PT% prior to lyophilization. Two lots of the plasma plus rFVII had a PT% of about 100% prior to lyophilization. After lyophilization, the PT% was about 85%. The PT% calibration curve from such reconstituted plasma was used to calculate PT% results. Values were very similar to those obtained using a calibration curve from Coag Cal N ("CCN") plasma, a lyophilized normal plasma containing all clotting factors.

Three lots of the calibrator plasma were produced by adding rFVIIa to a pool of HEPES buffered citrated plasma. The accelerated stability studies showed that after 35 days at 37°C (equivalent to 2 years at 4°C), the results were similar to those of CCN plasma and suggests they will have a similar stability. In two lots, the PT% was adjusted to approximately 100%

10

before lyophilization. Lyophilization appeared to reduce the PT% to between 90-95%. The pre-lyophilization target for the third lot was changed to between 105% and 108%, inclusive. Post-lyophilization, the third lot had a PT% of about 100%. The reconstituted third lot was stable for 8 hours at 4°C and room temperature. The PT% calibration curves from such lot were stable for 30 minutes.

As more fully explained in the examples, one or more of the following reagents were used in the examples that follow. (These examples are intended for purposes of illustration of the invention, not for limitation of the invention. For instance, the addition of HEPES is referred to as "dropwise" in an example. The invention obviously is not limited to use of the HEPES buffer or its dropwise addition.)

## Recombinant material:

Material	Lot No.	Concentration	Source
rFVIIa	29491	1 mg/ml (Novo)	Düdingen
rFVIIa	8293	1.2 mg/ml (Novo)	Harrow
rFVII	28193	30-40 U/ml	Harrow
rFVII	9393	24 U/ml	Harrow
rFVII	10393	15 U/ml	Harrow
rFVIIa	21593	2500 U/ml	Harrow

## 11

Harrow refers to the Haemostasis Research Group,  
Clinical Research Centre, Watford Road, Harrow,  
Middlesex, England.

## Other reagents:

- 5 TIS                      Thromboplastin IS lots TPS - 46 and 59  
                              (Baxter's dried rabbit brain with  
                              calcium PT assay reagent)
- 10 Innovin<sup>™</sup>                Innovin<sup>™</sup> reagent lots TFS - 12, 13, 14  
    and 24
- Saline                NaCl (0.9%) lots H1-75
- 15 Owrens Buffer    Owrens Buffer lots 550.029, 550.030 and  
    550.032
- Factor VII  
    Immuno Absorbed  
20 Plasma ("IAP") Factor VII IAP lots IAP7-25A and 26A
- FVII(a)  
    -Tris Buffer    Tris Buffer pH 7.4 lots H1-85  
                              Buffer used to dilute rFVII  
25                               (Although TRIS buffer is used in these  
                              examples, it is believed that any  
                              buffer of the same pH can be used.)  
                              0.05M Tris (hydroxymethyl)-aminomethane  
                              0.15M NaCl
- 30        Several lots of CCN plasma, a lyophilized normal  
          plasma containing all clotting factors, were tested  
          for PT% using TIS and Innovin<sup>™</sup> reagents. They were  
          also tested for the FVII% level. The results are  
          tabulated below. The five lots of CCN plasma were  
35 combined to make FNP 870.003.

12

	CCN lot No.	PT% TIS	PT% Innovin™	FVII% level
	540.042	92	—	98
	540.049	91	85	105
	540.050	100	85	—
5	540.053	97	85	97
	540.054	—	—	—
	FNP 870.003	100	100	100

Machines and software:

10

MLA Electra 1000C: No 572 - Software Version 3 Rev. E  
MLA Electra 900C: No 1753 - Software Version 4 Rev. 1  
MLA Electra 1000C: Software Version 5.0: Munchen

15

Methods:

## Prothrombin Time (PT)

The PT testing assays were performed as per the  
20 Box Inserts for TIS and Innovin™ reagents, and the MLA  
Electra 900C or 1000C operating manuals.

## Factor VII assay

Factor VII assays were performed as per the Box  
Insert of the FVII IAP and the MLA Electra 900C or  
25 1000C operating manuals. The dilutions of the plasma  
or concentrate were selected so that the clotting  
times obtained were within the range obtained using  
the calibration curve dilutions. In general, the 1 in  
10 dilution was assigned as 100% Factor  
30 VII.

EXAMPLE I

Recombinant FVIIa lot 29491 was diluted in Owrens Buffer  $1/10^2$ ,  $1/10^3$ ,  $1/10^4$ ,  $1/10^5$ ,  $1/10^6$ , and  $1/10^7$ .

Five 500 ul aliquots of FNP 870.003 were prepared. To 5 each of the aliquots of FNP 870.003 were added one 20 ul aliquot of one rFVIIa dilution. The PT% of the resulting plasmas were tested using the MLA Electra 900C.

When measured using TIS and Innovin<sup>™</sup> reagents, it 10 was possible to reduce the PT clotting time of FNP, thus increasing the % PT. (See Table 1a). Using the  $1/10^4$  dilution of the rFVIIa, the Factor VII% level in the FNP was raised by 13-20%. The calibration curves (Clotting Time, in seconds, vs.  $1/PT\%$ ) of FNP and FNP 15 plus rFVIIa were nearly parallel, indicating that the modified plasma (FNP plus rFVIIa) can be used as a calibrator. See Table 1b and Fig. 5.

Table 1: Effect of different concentrations rFVIIa in FNP on PT

20

Table 1a: rFVIIa dilutions in FNP

	PT Clotting Time (seconds)	
	TIS	Innovin <sup>™</sup>
25 None	14.4	11.4
10 2	9.8	8.5
10 3	10.1	9.7
10 4	13.6	11.1

14

Table 1a continued

10 5	14.6	11.5
10 6	14.7	11.6
Buffer	14.8	11.5

5

Table 1b: Data for rFVIIa 10 3 dilution in FNP, calibration curve, compared with data for FNP curve

Dilution	PT Clotting Time (seconds)			
	FNP		FNP + rFVIIa 10 3	
	TIS	Innovin™	TIS	Innovin™
Neat	14.5	11.5	12.6	10.5
1 in 2	20.2	15.0	17.4	13.7
1 in 4	32.0	25.5	27.9	22.4
1 in 8	56.1	44.3	51.7	41.2

10

15

EXAMPLE II

20 Recombinant FVIIa Lot 8293 was diluted in CCN plasma lot 049 by adding 50 ul of concentrated rFVIIa to 5 ml of CCN plasma, resulting in a 1 in 100 dilution, noted as 10 2). Then a range of 1 in 10 dilutions were produced by adding 500 ul of the

25 resulting plasma to 4.5 ml of the CCN plasma. Three further dilutions were made, resulting in  $1/10^3$ ,  $1/10^4$  and  $1/10^5$  dilutions. The CCN plasma lots and the four dilutions were tested using TIS and Innovin™ reagents. The results are set forth below:

30



Table 2: Addition of rFVIIa to CoagCal N plasma

Sample	TIS				Innovin™			
	Neat	1 in 2	1 in 4	1 in 8	Neat	1 in 2	1 in 4	1 in 8
CCN 042	15.1	22.5	40.2	73.4	11.9	16.3	27.3	53.1
CCN 049	15.0	22.6	37.6	75.6	12.0	16.0	27.5	51.1
10 5	12.2	17.3	29.7	57.6	10.3	13.5	21.8	40.6
10 4	10.7	15.4	25.6	50	9.7	12.2	19.1	35.6
10 3	9.8	13.6	22.0	41.5	8.9	11.0	16.3	28.5
10 2	9.5	13.2	21.9	42.6	—	10.9	16.3	—

CCN plasma, like FNP, experienced a reduction in PT, thus increasing the % PT by the addition of rFVIIa.

#### 15 EXAMPLE III

Testing was done on rFVII material. Reagents included FNP 870.003 and CCN plasma lot 042. Testing was performed on the MLA Electra 1000C.

Different volumes of the three lots of rFVII were added to CCN plasma lot 042. Because the rFVII preparation had lower FVII activity than the rFVIIa preparation, instead of diluting the FVII preparation and adding the dilution to the plasma as in Example I, a different method was used as described below. This was done by reducing the amount of distilled water added to reconstitute the CCN plasma by the volume of rFVII added. For example, when 100 ul rFVII was added, the vial of CCN plasma was reconstituted with only 900 ul of distilled water.

16

The % PT and FVII % of the reconstituted CCN plasma lot 042 samples were calculated using the FNP calibration curve assigned as 100% PT activity. The mean % PT for all calibration curve dilutions was used. All lots of rFVII raised the % PT and the Factor VII% levels. The effect on the % PT was not proportional to the rise in FVII activity. It was thought that the lots of "rFVII" may have activated rFVII present in variable amounts which led to a variable effect on the % PT which was not related to the assigned Factor VII level. Because of the variability in rFVII, the use of rFVIIa would be preferable as it would be a more consistent reagent.

Table 3: Addition of rFVII to CoagCal N plasma

15

Sample		PT%		Factor VII%	
		TIS	Innovin™	TIS	Innovin™
FNP CCN		100 90	100 95	— 98	— 98
rFVII 9393	75 ul	92	98	112	110
	100 ul	97	102	118	116
	100 ul	96	100	130	124
rFVII 28192	20 ul	108	108	116	104
	30 ul	108	115	200	202
	50 ul	121	120	284	304
rFVII 10393	150 ul	94	93	138	116
	200 ul	97	94	134	118

25

Factor VII calculated using CCN plasma as calibrator PT % calculated assuming FNP = 100%

30

EXAMPLE IV

The method of measuring the Factor VII in the concentrate was investigated and the relationship of the Factor VII levels and PT% in CCN plasma with different amounts of rFVIIa added was examined. Reagents included rFVIIa lot 21593, CCN plasma lots 042 and 049 and IAP7-26A.

Testing was performed on the MLA Electra 1000C. The Factor VII level of the rFVIIa preparation was measured in two ways, by adding the preparation to CCN plasma and assaying dilutions of 1/100 to 1/1000 in Owrens buffer.

A primary dilution of rFVIIa in CCN plasma was made (CCN plasma 5ml plus 20 ul rFVIIa), also referred to as "plasma + rFVIIa". Then the following dilutions were made from the plasma + rFVIIa and CCN plasma. See Table 4. The 10 ul dilution of Table 4, marked with the "+", is the same 10 ul dilution used in Table 5.

Table 4: Dilutions of rFVIIa in CoagCal N plasma

Amount of rFVIIa in 5 ml CCN	10 ul*	5 ul	4 ul	3 ul	2 ul	1 ul	0 ul
Primary dilution	1 ml	0.5 ml	0.4 ml	0.3 ml	0.2 ml	0.1 ml	0 ml
CCN	1 ml	1.5 ml	1.6 ml	1.7 ml	1.8 ml	1.9 ml	2.0 ml

Table 5: Further dilutions of rFVIIa

5	Amount of rFVIIa in 5ml CCN	2 ul	1 ul	0.5 ml	0.25 ul	0.125 ul	0 ul
	10 ul dilution*	1 ml	0 ml	0 ml	0 ml	0 ml	0 ml
10	CCN	4 ml	1 ml	1 ml	1 ml	1 ml	2 ml
	Mix + Transfer Previous Dilution	0 ml	1 ml	1 ml	1 ml	1 ml	0 ml

15

The results are found in Tables 8 and 9. The following formula was used to calculate the FVII% concentration in U/ml.

20

$$\frac{\text{FVII\%}}{100} \times 5 - 5 \times \left( \frac{1000}{\text{ul of rFVIIa added}} \right) = \text{FVII of the concentrate (U/ml)}$$

25

FVII%/100

100% FVII = 1 U/ml

x 5

5 ml of plasma

- 5

5 U/ml of FVII in this 5 ml of normal plasma

30

1000/ul

rFVIIa added

Volume of rFVIIa compared to 1000 ul added

35

Results were calculated for plasma + rFVIIa using the formula set forth above.

19

Table 6: Calculation of FVII levels - rFVIIa added to plasma

5	Amount rFVIIa added	rFVIIa FVII%	FVII (U/ml)
	5 $\mu$ l	253	1532
	4 $\mu$ l	246	1830
	3 $\mu$ l	224	2070
10	2 $\mu$ l	197	2437
	1 $\mu$ l	167	3345
	mean	218	2243

Table 6 shows that the concentrate FVII level was 2243 U/ml when rFVII was added to plasma.

15 Table 7: Calculation of FVII levels - dilutions of rFVIIa in buffer

	Dilution	FVII%	FVII% effective	FVII (U/ ml)
20	1/1000	441	44100	
	1/2000	292	58400	
	1/4000	190	76000	
	mean		59500	595
25	1/10000	93	93000	
	1/20000	51	102000	
	1/40000	29	116000	
	mean		103667	1037

Table 7 shows that the concentrate FVII level was between about 600 and 1000 U/ml when rFVII diluted in 30 buffer was tested.

When measuring rFVIIa in plasma, the result obtained (2243 U/ml) was similar to the quoted concentrated from Harrow (2500 U/ml). Estimates using diluted concentrate were lower (600-1000 U/ml) and we 35 concluded that this method is not useful.

A progressive rise occurs in % PT and FVII levels with increasing the addition volume of rFVIIa to plasma (Tables 8 and 9). As seen from the data in Table 9, there was a relationship between the rise in Factor VII level and the rise in PT%,  $r = 0.9661$ .

Table 8: Effect of rFVIIa on Prothrombin Time  
Lot rFVIIa 21593

10	Amount rFVIIa Added	Innovin™				Test Mode
		PT% Calibrator curve dilutions				
		Neat	1 in 2	1 in 4	1 in 8	
15	10 ul	10.1	14.8	24.2	50.4	10.5
	5 ul	10.4	15.3	25.5	49.7	10.5
	4 ul	10.6	15.7	26.7	53.8	10.9
	3 ul	10.8	16.5	28.3	55.9	10.8
	2 ul	11.1	16.4	28.2	58.5	11.3
	1 ul	11.5	17.4	30.8	60.6	11.4
	zero	12.4	19.2	34.2	67.2	12.5
20	2 ul	10.8	15.8	27.4	54.6	10.9
	1 ul	11.0	16.7	29.2	61.5	11.0
	0.5 ul	11.3	17.2	30.8	64.0	11.4
	0.25 ul	11.6	17.7	31.6	64.6	11.7
	0.125 ul	11.7	18.1	32.8	66.0	11.7
	zero	12.0	18.9	34.4	62.6	12.0

	Amount rFVIIa Added	Thromboplastin IS				Test Mode
		PT% Calibrator curve dilutions				
		Neat	1 in 2	1 in 4	1 in 8	
30	10 ul	11.4	17.7	30.5	62.1	—
	5 ul	12.0	19.0	32.0	66.9	—
	4 ul	12.3	19.0	33.6	66.5	—
	3 ul	12.4	19.3	33.7	71.4	—
	2 ul	12.8	20.4	34.4	73.0	—
	1 ul	13.3	21.2	36.3	76.3	—
35	zero	15.2	24.7	41.9	85.9	—

21

Table 9: Investigation of increasing PT% and FVII% level

Innovin™				
Amount rFVIIa Added	FVII% Dilution 1 in 10	PT%	Rise in PT%	Rise in FVII%
10 ul	—	131	33	—
5 ul	253	131	33	142
4 ul	246	123	25	135
3 ul	224	125	27	113
2 ul	197	115	17	86
1 ul	167	114	16	56
zero	111	98	0	0
2 ul	—	118	18	—
1 ul	—	116	16	—
0.5 ul	—	109	9	—
0.25 ul	—	104	4	—
0.125 ul	—	104	4	—
zero	—	100	0	—

Table 9 continued

Thromboplastin IS				
Amount rFVIIa Added	FVII% Dilution 1 in 10	PT%	Rise in PT%	Rise in FVII%
10 ul	—	161	61	—
5 ul	—	146	46	—
4 ul	—	140	40	—
3 ul	—	138	38	—
2 ul	—	131	31	—
1 ul	—	123	23	—
zero	—	99.5	0	—

EXAMPLE V

Stable lyophilized plasma which has had Factor rFVII added to be used as a calibrator in the PT % test was prepared as follows.

Reagents used were rFVII lots 28193 and 9393, plasma as described in Table 11a, and HEPES buffer H1-83.

The volume of rFVII needed was calculated as follows. It was expected that after lyophilization the PT % would be about 85-90%; thus a rise in PT % of 10-15% was required. Preliminary work with lot 28193 suggested that 20-30 ul rFVII per 5 ml of plasma created the desired rise in PT %; 25 ul rFVII per 5 ml of plasma was used. Lot 9393 had a lower Factor VII level and about 400 ul rFVII per 5 ml plasma was needed to raise the PT%.

Ten units of approximately 200 ml each of plasma were selected from each of the plasma bags described in Table 11a. All plasma had been collected into the anticoagulant CPD-A. The plasma was carefully thawed in a large waterbath at 37°C. Bag contents were mixed until all ice had disappeared and the plasma was free from undissolved precipitate. Once thawed, the bags were kept in crushed ice. The entire contents of each bag were pooled and stirred thoroughly while kept cool by crushed ice. Four pools were prepared for lyophilization as described in Table 10.



Table 10: Preparation of different plasma pools

Pool	Volume of plasma	Volume of HEPES	Volume recombinant FVII
Pool P	100 ml	None	None
Pool P1	100 ml	3 ml	None
Pool P1 & 28193	100 ml	3 ml	500 ul of lot 28193
Pool P1 & 9393	100 ml	3ml	8 ml of lot 9393

HEPES was added dropwise to the stirred plasma.

10 The recombinant FVII was added last and the final mixture stirred thoroughly. The resulting plasmas were pipetted into separate 1.1 ml vials and stored at 4-8°C for about 1 hour before lyophilization. After lyophilization, the vials were kept at 4-8°C. Several 15 vials from each lot were not lyophilized but stored at -70°C storage. Prior to lyophilization, the different pools of plasma, CCN plasma lot 042 and CCN plasma lot 049, both freshly reconstituted, were tested using the MLA Electra 1000C. Both PT% and FVII% level assays 20 were performed. Plasma samples were tested in the calibration curve mode and in the test mode. The 10 plasmas used to make up Pool P were all normal (see Table 11a). Testing of fresh pools suggest a PT% of approximately 100% in both pools (P1 & 28193, P1 & 25 9393). The rise in FVII was 140% and 190%, respectively. The conclusion is that the amount of rFVII needed to prepare a control with approximately

24

100% PT can be predicted. Prior to lyophilization, the addition of HEPES buffer reduced the PT% by about 5%.

After lyophilization, the different pools of plasma, CCN plasma lot 042 and CCN plasma lot 049 were tested using the same instrument and procedure as in their testing before lyophilization. After lyophilization, the pool without HEPES showed a loss of 11% PT whereas the pool with HEPES showed no difference. The two pools with rFVII showed a slight (2%) loss in % PT. No changes were seen in FVII% levels after lyophilization even in the pool without HEPES. (See Tables 12a-12b.)

Table 11a: PT Clotting Time of plasmas making up plasma Pool P

BAG No.	Clotting Time	
	1020 hrs	1344 hrs
1	11.6	11.2
2	11.5	11.1
3	12.8	12.4
4	11.8	11.4
5	12.7	12.3
6	11.8	11.3
7	12.3	11.6
8	12.3	11.7
9	12.1	11.6
10	11.2	11.0
Mean	12.01	11.56

25

Table 11b: PT Clotting Times before  
lyophilization Innovin™ PT Reagent

Calibrator	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
CCN 049	12.2	18.0	32.7	69.0	12.14	84.2
CCN 042	12.2	18.3	31.9	68.3	12.16	84.0
Pool P	11.6	17.6	32	66.9	11.56	91.3
Pool P1	12.1	18.6	32.8	66.2	11.99	86.0
Pool P1 & 28	11.1	15.7	28.4	58.1	10.95	100.1
Pool P1 & 93	11.0	16.1	27.0	58.6	11.06	98.4

Table 11c: Factor VII assay before lyophilization  
Innovin™ PT Reagent

Calibrator	Calibration Curve					Test Mode 1 in 10	FVII% 1/10**	FVII% 1/20**
	1 in 10	1 in 20	1 in 40	1 in 80	1 in 160			
CCN 042	23.1	21.4	41.2	53.6	69.7	23.2	110	104
CCN 049	23.4	31.5	42.1	54.1	70.4	23.1	111	103
Pool P	24.3	32.9	43.9	57.7	74.3	23.5	107	93
Pool P1	24.5	33.9	43.6	56.5	72.4	23.6	105	86
Pool P1 & 28	20.4	27.8	35.8	47.4	63.3	20.8	144	141
Pool P1 & 93	18.1	23.6	NA	39.8	54.1	18.3	198	210

\*\* : Calculated with CCN plasma 049 calibration Curve:

PT% = 85%, FVII% = 105%

"Test Mode" means that a sample can be tested as a calibrator whether the sample is diluted or just as neat plasma.

Table 12a: PT Clotting time after lyophilization  
Innovin™ PT Reagent

Innovin™	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
CCN 049	12	17.5	30.8	63.1	12.15	84.1
CCN 042	12.2	17.8	30.9	65.7	12.05	85.3
Pool P	12.5	18.9	33.1	66.9	12.55	79.9
Pool P1	11.8	18.1	31	65.4	11.85	87.6
Pool P1 & 28	11	16.3	27.4	57.9	11.05	98.5
Pool P1 & 93	11.1	15.9	26.5	54.6	11.2	96.3

26

Table 12b: Factor VII assay after lyophilization

	Calibrator	Calibration Curve					Test Mode		FVII% **
		1 ln 10	1 ln 20	1 ln 40	1 ln 80	1 ln 160	1 ln 10	1 ln 10	
5	CCN 049	23.1	30.2	40.9	54.2	71.7	23.5	23	112
	CCN 049	23.5	31.1	41.1	53.6	71.3	—	—	—
	CCN 042	23.6	31.7	42.9	56.1	72.6	23.6	23.6	105
	Pool P	—	—	—	—	—	23.2	24.1	100
	Pool P1	—	—	—	—	—	23.9	23	112
10	Pool P1 & 28193	—	—	—	—	—	20.5	20.3	153
	Pool P1 & 9393	—	—	—	—	—	18.2	18.1	203

\*\* : Calculated with CCN plasma 049 calibration Curve:  
PT% = 85%, FVII% = 105%

#### EXAMPLE VI

Stable, lyophilized plasma to which rFVIIa has been added to be used as a calibrator in the PT % test was prepared as follows.

Reagents used were rFVIIa Lot 21593, Pool 2 (CCN plasma lot 053 just before lyophilization), and TRIS Buffer Lot H1-85. Four hundred milliliters of a plasma pool ready to use (containing HEPES) were used to prepare CCN plasma lot 053.

Table 13: Preparation of different plasma pools

Pool Name	rFVIIa added	Plasma pool 053	[FVII] added
Pool P2	None	100 ml	None
Pool P2/20	20 ul	100 ml	1 ul/5ml
Pool P2/10	5 ul	50 ml	0.5 ul/5ml

Pool Name	Tris Buffer added	Plasma pool 053	[Tris B.] added
Pool P2/B	20 ul	100 ml	1 ul/5ml

27

Vials were filled with 1.1 ml pooled plasma and stored at -70°C for five days and then lyophilized. A certain number of vials were kept at 70°C and not lyophilized.

5 Prior to lyophilization, the four pools were tested on the MLA Electra 1000C. After lyophilization, the four pools were tested against the corresponding four frozen pools using the same instrument and procedure as in the testing prior to  
10 lyophilization.

The results found in Tables 14 and 15 were calculated using a previous CCN plasma lot 049 calibration curve (Table 11b: 12.2, 18.0, 32.7 and 69.0 seconds assigned as 85% PT). The fresh results  
15 for two lots (Pool P2/20 and Pool P2/10) are 102% and 98% respectively. After lyophilization, there appears to be a 5-8% drop in the PT %, which did not occur with rFVII. The frozen samples did not show this drop. It appears that in plasmas where rFVIIa is used  
20 to increase the PT%, there is a 5-8% loss of PT% during lyophilization. This needs to be compensated for during the manufacturing process.

Table 14: Pool P2 fresh before the lyophilization

Innovin™	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
Pool P2	12.3	17.7	30.5	60.1	12.1	84.7
Pool P2/B	12.1	17.6	30.1	59	12	85.5
Pool P2/10	11.1	15.6	26	51.1	—	97.8
Pool P2/20	10.8	15.2	24.5	50.2	10.8	102.5

10 \*\* Calculated with CCN plasma lot 049 Calibration Curve: PT% = 85%, FVII% = 105%

Table 15: Pool P2 after the lyophilization

Lyophilized calibrator test with Innovin™ PT Reagent						
Calibrator	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
Pool P2	12.2	18.6	32.3	66.6	12.3	82.5
Pool P2/B	12	18.6	32.1	65.1	12.3	82.5
Pool P2/10	11.4	16.7	29.1	60.3	11.65	89.5
Pool P2/20	11.1	16.6	28.6	57.7	11.3	94.8
Frozen calibrator test with Innovin™ PT Reagent						
Calibrator	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
Pool P2	11.7	17.4	30.2	59.3	11.65	90.1
Pool P2/B	11.7	17.2	30.1	59.4	11.85	87.6
Pool P2/10	11	15.5	26.3	52	11.1	97.8
Pool P2/20	10.7	15.3	25.5	51	10.5	107.7
Lyophilized calibrator test with TIS PT Reagent						
Calibrator	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
Pool P2	15.2	20.1	34.3	66.2	15.7	—
Pool P2/B	14.9	22.7	37.7	73.3	15	—
Pool P2/10	14.1	21.1	34.3	67.8	14.1	—
Pool P2/20	13.6	20.1	34.3	66.2	13.7	—

Table 16: Frozen Pool P2

Frozen calibrator test with TIS PT Reagent						
Calibrator	Calibration Curve				Test Mode	PT% **
	Neat	1 in 2	1 in 4	1 in 8		
5 Pool P2	14.7	22.6	37.5	74.3	14.9	—
Pool P2/B	14.9	22.5	37.3	72.6	15	—
Pool P2/10	13.7	20.3	33.8	67.6	13.7	—
Pool P2/20	13.1	19.5	32.9	64.6	13.1	—

10 \*\*: Calculated with CCN plasma lot 049 calibration curve: PT% = 85%, FVII% = 105%

#### EXAMPLE VII:

15       The accelerated stability of Pool P2 (as prepared in Example VI) with rFVIIa added was tested and compared with two lots of CCN plasma. Reagents used were Pools P2, P2/13, P2/10, P2/20, CNN plasma lots 050 and 053. Several vials of the plasmas were stored

20 at 37°C and tested after 10, 14, 26 and 35 days on the MLA Electra 1000C according to the Box Insert and the MLA Electra 1000C Handbook. Vials of the same plasmas stored at 4°C were tested for the same time periods. All plasma tested showed a progressive drop in the PT%

25 on incubation at 37°C. The plasma containing rFVIIa did not drop differently than those not containing rFVIIa. Adding rFVIIa does not change the stability of the plasma incubated at 37°C measured using the PT% assay. See Table 17. Subsequent analysis of further

30 lots with Arrhenius stability testing has given a predicted shelf life of greater than 2 years.

Table 17: Accelerated stability Pool P2

5	Callibrator	Prothrombin Time in %							
		10 days		14 days		26 days		35 days	
		4°C	37°C	4°C	37°C	4°C	37°C	4°C	37°C
	CCN 050	83.6	79.4	84.7	76.5	83.6	76.5	83.6	79.7
	CCN 053	83.6	78.4	83.6	76.5	82.5	75.7	83.6	73.0
	Pool P2	83.6	78.4	83.6	76.5	83.6	73.8	83.6	71.3
	Pool P2/B	84.7	78.4	84.7	76.5	83.6	73.8	84.7	73.0
10	Pool P2/10	92.1	85.8	90.8	83.6	92.1	80.4	92.1	79.4
	Pool P2/20	94.8	89.5	96.3	85.8	94.8	83.6	96.3	81.4

PT% is calculated with CCN lot 049 Calibration Curve:  
PT% = 85%.

15

EXAMPLE VIII:

A previously prepared pool of citrated plasma, from 10 donors (See Table 6a), stored at -20°C, was thawed in a 37°C waterbath and then stored at 4°C.

20 When the temperature of the thawed plasma reached 4°C, then a HEPES solution was added slowly dropwise.

The HEPES solution was prepared by adding 40 mg of HEPES powder to 100 ml distilled water. The pH was adjusted to approximately 7.3 to 7.5 using 5M NaOH.

25 (About 5 ml of 5 M NaOH was needed.) This resulted in a 40% HEPES solution (lot H1-83).

For each liter of plasma in the pool, 30 ml of the 40% HEPES solution were added. The pooled plasma and the HEPES solution were mixed for 10 minutes, with  
30 care not to create foam.



31

All testing was performed using an MLA Electra 1000C. The PT% of the pool plus HEPES buffer (the "Buffered Pool") was determined. Recombinant Factor VIIa was then added to the Buffered Pool in a step-wise manner, as described below, until the PT% of the Buffered Pool plus rFVIIa was between 105% and 108%. It was adjusted 5-8% above 100% to allow for PT% loss of 5-8% during lyophilization. The rFVIIa had previously had its activity determined by adding dilutions to plasma (as described in Example IV), and this activity was used in the following formula to determine the amount (in ml) of rFVII to add per ml of Buffered Pool.

Amount rFVIIa (ml) = 
$$\frac{\text{Volume of Buffered Pool (ml)} \times \text{U/ml rFVII required}}{\text{rFVIIa concentration (U/ml)}}$$

The total amount of rFVIIa that should be added to the Buffered Pool to achieve a PT% of 105-108% is about 0.6 Units per ml of plasma. If the PT% is as follows, then the amount of rFVIIa that is required is as follows:

< 90%	add 0.6 U/ml
< 100%	add 0.3 U/ml
< 105%	add 0.15 U/ml.

Once the target PT% activity of the plasma pool with rFVIIa was achieved, the mixture was again thoroughly stirred for at least two minutes. Two

32

aliquots of the mixture were tested and the mean of all 8 results was calculated. If the mean result was between 105-108% (inclusive), the material was accepted for lyophilization. If need be, further buffered plasma that has not had rFVIIa added to it can be added to the Buffered Pool to reduce the PT% to achieve the required value. See Tables 18-20. Data from a pilot production size run is shown in Tables 18, 19, and 20.

10            Table 18: Pre-lyophilization Testing  
Calibration plasma curve

	Dilution of Plasma			
	Neat	1/2	1/4	1/8
15      Sample 1	11.8	16.8	27.9	56.7
	11.8	16.6	27.7	54.5
Sample 2	11.5	16.7	29.8	55.7
	11.3	16.2	28.3	56.5
20      Sample 3	11.4	16.3	27.2	54.6
	11.3	16.1	27.7	53.9
Mean*	11.5	16.5	28.1	55.3
PT%	88	44	22	11
			COD	0.939

25

Table 19: Reagents used

5	Calibration Plasma	
	Lot No.	R & D Pool P3
	PT% with Innovin	88
10	rFVIIa conc. Lot No.	21593
	rFVIIa conc. (U/ml)	2500
	Innovin™ Lot No.	TFS-12
15	Saline Lot No.	H1-86
	Machine type	1000C
	Machine No.	187
	Programme version	3.E

Testing of Plasma Pool:

20	Measured volume (ml)	<u>1200</u>
	Reserved plasma volumes (ml)	<u>200</u>
25	Numbers of donor units	<u>10</u>

Table 20: Results obtained

	Material Tested	Raw Data			
		Clotting Time (Secs)		Calculated PT%	Mean
30	Initial Pool R&D Pool P3 Calibration Plasma	12.5	12.3	75.6	77.8
		12.2	12.5	78.9	75.6
		12.3	12.4	77.8	76.7
35	Plasma Pool plus 0.3 U/ml plasma rFVIIa Volume of rFVIIa added = 0.12 ml	11.2	10.9	92.5	97.5
		11.0	11.1	95.8	94.1

Table 20 Continued

5	Further addition of							
	rFVIIa		Plasma (ml)					
	Volume U/ml Plasma							
	0.06	0.15	0	95.8	99.3	11.0	10.8	98.9
				99.3	101.2	10.8	10.7	
	0.06	0.15	0	101.2	105.1	10.7	10.5	104.7
				105.1	107.2	10.5	10.4	
	0.03	0.075	0	105.1	105.1	10.5	10.5	104.7
				101.2	107.2	10.7	10.4	
	0.03	0.075	0	107.2	107.2	10.4	10.4	104.8
107.2				109.4	10.4	10.3		
10	0	0	103.2	109.4	10.6	10.3	107.3	
			107.3	109.4	10.4	10.3		
Final PT%			107.3					

Testing of the lyophilized product was performed using the MLA Electra 1000C. A lyophilized plasma that had been calibrated against FNP was used as a calibrator ("the Calibrator"). The PT% of the lyophilization product was calculated (using only the results from the undiluted plasma). This was assigned as the PT% of the product. A calibration curve was then obtained using the lyophilized product. As an in-process control check of the lyophilized product, a range of test results (see Tables 21-24) were calculated using the lyophilized product and the

Calibrator, and the percentage difference was calculated. The lyophilized product was deemed to be acceptable if there were no differences greater than 15% (See Table 25).

5           **Table 21: Post-Lyophilization Testing:  
Calibration plasma curve**

		Dilution of Plasma			
		Neat	1/2	1/4	1/8
10	Sample 1	11.5	16.3	27.9	53.5
		11.5	16.2	26.4	53.7
	Sample 2	11.3	16.1	27.8	53.0
		11.2	16.0	26.7	54.6
	Sample 3	11.2	16.1	27.5	52.8
		11.2	16.1	26.6	53.3
	Mean	11.5	16.1	27.2	53.5
	PT%	86	44	22	11
15				COD	0.931

**Table 22: Reagents Used**

20	Calibration Plasma	
	Lot No.	R&D Pool P3
	PT% with Innovin™	88
	Innovin™ Lot No.	TFS-12
25	Saline Lot No.	H1-86
30	Machine Type	1000C
	Machine No.	187
	Programme version	3.E

Table 23: Innovin PT Calibrator:

	Dilution of Plasma			
	Neat	1/2	1/4	1/8
Sample 1	10.4	15.1	26.0	54.1
	10.6	15.1	26.5	51.2
Sample 2	10.6	15.5	26.8	51.1
	10.7	15.3	26.1	50.8
Sample 3	10.5	15.8	27.1	51.1
	10.7	15.5	26.8	52.1
Mean	10.6	15.4	26.5	51.2
PT%	104	52	26	13
			COD	0.947

Table 24: Final Results of Pool P4

IPTC Lot PT% with Innovin™	R&D Pool P4
	104
rFVIIa conc Lot No. Units added/ml plasma	21593 0.75

Table 25: Comparison of calculation with IPTC and CoagCal N:

Test Results (secs.)	PT% Calculated Using IPTC	PT% Calculated Using CoagCal N	% Difference CCPT/CCN
9	142.5	151.0	5.6
9.5	127.4	131.7	3.2
10	115.2	116.8	1.4
11	96.7	95.3	1.5
12	83.3	80.4	3.6
13	73.1	69.6	5.0
14	65.1	61.33	6.19

Table 25 Continued

5	16	53.6	49.6	8.1
	20	39.5	35.8	10.3
	25	29.7	26.6	11.7
	30	23.8	21.1	12.8
	40	17.1	15.0	14.0
10	50	13.2	11.6	13.8
	60	10.9	9.5	14.7
	70	9.2	8.0	15.0
	Final PT%			104

EXAMPLE IX

The stability of the first Pilot production  
 15 (identified here as Lot P4), when reconstituted, was  
 tested using Innovin<sup>™</sup> reagent lot TFS-12 and Saline  
 H1-86. Dilution stability testing was performed on  
 the MLA Electra 900C using programme version 4.1.

The testing was performed immediately after the  
 20 dilutions had been prepared ( $t = 0$ ) and exactly 30  
 minutes after preparation ( $t = 30$ ). All reconstituted  
 stability testing was performed using the E1000C using  
 programme version 3E. Six vials were reconstituted.  
 Three were left at room temperature for 8 hours and  
 25 three at 4°C for 8 hours. After 8 hours, three more  
 vials were freshly reconstituted and all nine vials  
 had calibration curves produced. Innovin<sup>™</sup> reagent was  
 freshly reconstituted and tested immediately, after 4

and 8 hours stored on the E1000C (8°C) using freshly reconstituted reagents at each time point.

The material was deemed not to have failed stability testing if the clotting time in seconds was not more than 10% different from the clotting time obtained from lyophilized material stored at 4°C that had been freshly tested after reconstitution. The plasma dilutions were stable for 30 minutes. See Table 26. There was no significant variation between time 0 and time 30. The stability of Innovin™ reagent on the Electra 1000C is also good for 8 hours at 8°C (Table 27). The results do not give a variation from time  $t = 0$  until time  $t = 8$  hours. The PT calibrator reconstituted stability was also measured for 8 hours. There is little change between time  $t = 0$  and time  $t = 8$  hours at either 4°C (2%) or room temperature (4%). See Table 28.

The testing confirms the stability of the dilutions, the stability of Innovin™ TFS-12 and the reconstituted stability of Pilot lot P4.



**Table 26: Stability of the Dilutions of PT Calibrator:**

	Incubation Time							
	0 Minutes				30 Minutes			
	Dilution of plasma				Dilution of Plasma			
	Neat	1 in 2	1 in 4	1 in 8	Neat	1 in 2	1 in 4	1 in 8
Sample 1	11.1	13.5	20.0	36.0	11.1	13.1	20.9	37.0
	10.9	13.0	19.2	36.6	10.8	13.0	19.3	35.9
Sample 2	11.4	13.0	19.4	37.9	11.0	12.7	19.7	37.4
	10.7	12.6	19.6	36.0	10.5	12.7	19.8	34.8
Sample 3	10.7	12.5	21.0	37.5	10.9	12.8	20.0	35.0
	10.4	12.5	19.1	38.4	10.4	12.9	18.9	24.5
Mean	10.87	12.85	19.72	37.07	10.78	12.87	19.77	35.77

Table 27: Stability of Innovin™ Reagent  
on the Electra 1000C

Sample	Incubation Time											
	0 Hours				4 Hours				8 hours			
	Dilution of plasma				Dilution of plasma				Dilution of plasma			
	Neat	1 In 2	1 In 4	1 In 8	Neat	1 In 2	1 In 4	1 In 8	Neat	1 In 2	1 In 4	1 In 8
1	10.9	15.3	24.9	49.0	10.7	15.0	25.3	49.9	10.8	16.0	27.3	53.9
	11.0	14.9	24.4	48.9	10.7	15.1	25.2	49.4	10.6	15.9	26.8	52.6
2	11.0	16.2	27.5	54.8	10.9	16.0	27.3	54.5	10.8	16.5	27.8	53.2
	10.7	16.2	26.9	55.2	11.0	15.9	27.8	55.5	10.7	16.0	26.7	53.4
3	10.8	16.0	28.1	52.7	10.9	16.2	27.6	55.5	10.7	16.6	27.3	54.0
	10.8	16.2	27.9	53.0	10.8	16.1	27.4	61.7	10.9	16.0	26.9	55.5
Mean	10.87	15.80	26.62	52.27	10.83	15.72	26.77	54.42	10.75	16.17	27.13	53.77

5

10

Table 28: PT Calibrator: reconstituted stability of 8 hours

Sample	Incubation Time											
	0 hours				4 hours				8 hours			
	Dilution of plasma				Dilution of plasma				Dilution of plasma			
	Neat	1 In 2	1 In 4	1 In 8	Neat	1 In 2	1 In 4	1 In 8	Neat	1 In 2	1 In 4	1 In 8
1	10.90	15.90	28.80	54.80	11.30	17.30	28.80	56.10	10.90	15.90	27.20	56.90
	10.70	16.00	27.30	56.00	11.20	16.50	28.60	55.90	10.80	15.90	27.50	52.30
2	10.90	15.90	28.10	54.20	11.30	16.10	28.70	NCD	11.10	16.10	27.80	54.40
	10.70	15.80	27.80	56.10	10.90	17.30	27.80	56.80	10.90	16.50	26.60	55.40
3	10.70	16.60	26.60	NCD	11.10	16.60	28.60	54.50	11.00	16.10	27.20	55.60
	10.50	16.00	27.50	53.50	10.90	16.50	27.40	55.10	10.80	15.70	28.60	52.90
Mean	10.73	16.03	27.35	54.92	11.12	16.72	28.32	55.68	10.92	16.03	27.52	54.58

RT = Room Temperature = 24°C. NCD = No Clot Detected

EXAMPLE X

Further stability testing was performed on Lot P4. The failure criterion was defined as a change of 10% in the Clotting Time (in seconds) as compared with the mean baseline value.

The accelerated stability calculation with the Arrhenius method was calculated with the SigmaPlot program as follows:

- 10        1. For each temperature plot decimal log of concentration (in this case - Clotting Time in seconds) (Y axis) against the time (in this case - days) (X axis) (Table 29).
2. For each temperature (graph) calculate the  
15 regression equation  $Y = m X + b$ .
3. Define a percent change at which the product is no longer acceptable, (in this case - +10% of Clotting time; mean baseline + 10% =  $10.69 + 10\% = 11.76$  seconds), convert the value of the zero time  
20 analyses to decimal log concentration (in this case - Log of Clotting Time(s) =  $\log$  of 11.76 - 1.070).
4. Using the regression equations for each temperature, substitute the decimal log and calculate the day failure.

5. Plot decimal log days from section 4, against 1/absolute temperature (Table 30).

6. Calculate the regression equation  $Y = m X + b$ , for the graph in section 5.

7. Using the regression equation from section 6, calculate the expected shelf life at 4°C.

Table 31 shows baseline data which demonstrates the reproducibility between different vials of Lot P4. Tables 32 and 33 show the results of testing of controls during the stability testing. Table 34 shows that the stability of Lot P4 failed after 45 days at room temperature (25°C). Table 35 shows the stability of Lot P4 at 30°C; Table 36 shows the stability of Lot P4 at 37°C; and Table 37 shows the stability of Lot P4 at 50°C.

Table 29: Calculation of failure day for each temperature

Temperature	X axis	Y axis		Statistics	Failure day
	Days	Clotting time (secs.)	Decimal Log of CT (secs.)		
at 25°C	0	10.69	1.029		
(Room temp.)	5	11.10	1.045		
	11	11.28	1.052		
	15	11.38	1.056		
	20	11.20	1.049		
	32	11.62	1.065		

Table 29 Continued

		37	11.52	1.061	$r = 0.933$	
		45	11.73	1.069	$l = 1.038$	
		56	12.00	1.079	$s = 0.00072$	44.4 days
5	at 30°C	0	10.69	1.029		
		3	11.13	1.046		
		4	11.47	1.060		
		7	11.65	1.066		
		8	11.63	1.066		
10		9	11.73	1.069		
		10	11.73	1.069		
		14	12.05	1.081		
		16	12.15	1.085	$r = 0.96956$	
		18	12.33	1.091	$l = 1.039$	
15		20	12.60	1.100	$s = 0.00304$	10.2 days
	at 37°C	0	10.69	1.029		
		1	11.30	1.053		
		2	11.45	1.059		
		3	11.62	1.065		
20		4	11.68	1.067	$r = 0.96541$	
		5	12.08	1.082	$l = 1.037$	
		6	12.32	1.091	$s = 0.009$	3.67 days

45

Table 29 Continued

5	at 50°C	0	10.69	1.029		
		0.083 (2 hrs)	11.12	1.046		
		0.167 (4 hrs)	11.32	1.054	r =	0.97223
		0.25 (6 hrs)	11.45	1.059	l =	1.033
		0.333 (8 hrs)	11.67	1.067	s =	0.1068
						0.346 day

Formula	Y = 1.070
Y = m X + b	m = slope (s)
X = Y - b/m	b = Intercept (l)

Table 30: Calculation of shelf life stability at 4°C

15	X-axis		Y axis		Statistics
	Temperature	1/temperature	failure day	log failure day	
	25°C	0.04	44.4	1.65	
	30°C	0.033	10.2	1.0086	r = 0.953
	37°C	0.027	3.67	0.56	l = 0.7046
20	50°C	0.02	0.346 (8.1 hours)	-0.46	s = 18.875

Formula	X = 0.25 (1.4°C)
Y = m X + b	m = Slope (s) = 18.875
	b = Intercept (l) = - 0.7046

46

$Y = 18.875 \times 0.25 + (-0.7046) = 4.014$   
 $\Rightarrow$  Inv. log of 4.014 = 10327 days stable  
 $\Rightarrow$  28.3 years  
 28 years - 33% = 18 years

5

In conclusion, Lot P4 is stable 18 years at 4°C

Table 31: Stability testing - Baseline Data

10

Assay: Prothrombin Time

Reagents: Innovin PT Calibrator lot  
(P4)

PILOT LOT 1

15

Innovin reagent lot

TFS - 12

Saline 0.9% Lot

8 H1 - 86

20

Machine: MLA E1000C  
Software Version

5.00E P46

Vials	Reference Tested	Clotting Time (Secs.)								
		Neat		1 in 2		1 in 4		1 in 8		
25	1	Fig. 14 NB CO82 P84	10.4	10.6	15.1	15.1	26.0	26.5	54.1	51.2
	2		10.6	10.7	15.5	15.3	26.8	26.1	51.1	50.8
	3		10.5	10.7	15.8	15.5	27.1	26.8	51.1	52.1
	4	Table 32 NB Co82 P 84	10.6	10.4	15.4	15.4	26.8	26.2	53.0	51.6
	5		10.6	10.5	15.5	15.4	26.1	26.1	52.3	52.6
	6		10.6	10.5	15.4	15.4	26.3	25.8	53.0	51.8
30	7	Table 45 NB Co95 p2	10.9	11.0	15.3	14.9	24.9	24.9	49.0	48.9
	8		11.0	10.7	16.2	16.2	27.5	26.9	54.8	55.2
	9		10.8	10.8	16.0	16.2	28.1	27.9	52.7	53.0



Table 31 Continued

10	Table 46 NB CO95 p3	10.9 10.7	15.9 16.0	26.8 27.3	54.8 56.0
11		10.9 10.7	15.9 15.8	28.1 27.8	54.2 56.1
5	Mean	10.69	15.60	26.67	52.70
	SD	0.178	0.382	0.902	1.996
	CV	1.666	2.447	3.384	3.787
10	Mean + 10%	11.76	17.16	29.34	57.97
	Mean - 10%	9.62	14.04	24.00	47.43

Table 32: Stability testing - Controls

Assay: Prothrombin Time

Reagents: Innovin™ lot  
Saline 0.9% lotTFS - 12H1-86 H1-87Machine: MLA E1000C  
Software version5.00 E P 46

25	Clotting Time (Seconds)				
	CoagCalN	Neat	1 in 2	1 in 4	1 in 8
	Sample 1	12.1	17.5	31.2	61.7
		12.1	17.5	29.5	56.1
	Sample 2	12.2	17.6	29.9	57.7
		12.0	17.9	30.0	58.6
	Sample 3	12.1	17.6	30.0	58.4
		12.0	17.7	30.1	57.6
30	Mean	12.08	17.63	30.12	58.35
	PT%	85	42.5	21.25	10.625

Coag Cal N

Lot No: 540.053

Innovin™ PT% 85

Table 32 Continued

	Date	Control	Clotting Time (seconds)		Mean	PT%
5	30.11.93	CTN	12.9	12.8	12.85	75.6
		CTP	21.0	20.8	20.90	36.0
	6.12.93	CTN	12.7	12.5	12.60	78.3
		CTP	20.7	21.5	21.10	35.5
10	17.12.93	CTN	12.9	12.7	12.80	76.1
		CTP	20.9	20.3	20.60	36.7
	20.12.93	CTN	13.1	12.9	13.00	74.1
		CTP	21.3	21.2	21.25	35.2
15	21.12.93	CTN	12.6	12.7	12.65	77.7
		CTP	20.5	20.4	20.40	37.2
	22.12.93	CTN	12.4	12.4	12.40	80.6
		CTP	20.3	20.4	20.35	37.3
20	23.12.93	CTN	12.6	12.4	12.35	81.2
		CTP	21.0	20.6	20.60	36.7
	27.12.93	CTN	12.6	12.5	12.50	79.4
		CTP	21.1	20.8	20.80	36.2
25	28.12.93	CTN	12.6	12.5	12.50	79.4
		CTP	20.8	20.5	20.50	37.0
30	29.12.93	CTN	12.5	12.5	12.50	79.4
		CTP	20.9	20.9	20.90	36.0

	Controls		Lot No.	Assigned Value
35	CoagTrol N	CTN	537.001	73 - 99%
	CoagTrol P	CTP	541.034	29 - 39%

49

Table 33: Stability testing - Controls

Assay: Prothrombin Time

5 Reagents: Innovin™ lot TFS - 12  
 Saline 0.9% lot H1-86 H1-87

Machine: MLA E100C Software version 5.00 E P 46

10

Clotting Time (Seconds)				
CoagCalN	Neat	1 in 2	1 in 4	1 in 8
Sample 1	12.1	17.5	31.2	61.7
	12.1	17.5	29.5	56.1
Sample 2	12.2	17.6	29.9	57.7
	12.0	17.9	30.0	58.6
Sample 3	12.1	17.6	30.0	58.4
	12.0	17.7	30.1	57.6
Mean	12.08	17.63	30.12	58.35
PT%	85	42.5	21.25	10.625

CoagCal N  
 Lot No: 540.053  
 Innovin™ PT% 85

15

Table 33 continued

20

Date	Control	Clotting Time(s)	Mean	PT%
30.12	CTN	12.5 12.3	12.40	80.6
	CTP	20.4 21.1	20.75	36.4
25 3.01.94	CTN	12.6 12.6	12.60	78.3
	CTP	21.3 22.3	21.80	34.0
5.01.94	CTN	12.6 12.5	12.55	78.8
	CTP	20.6 20.1	20.35	37.3
30 6.01.94	CTN	12.7 12.5	12.60	78.3
	CTP	20.7 21.1	20.90	36.0
06.01.94 MLA 900	CTN	12.5 12.5	12.50	79.4
	CTP	20.0 20.1	20.05	38.0

50

Table 33 Continued

5	7.01.94	CTN	12.7 12.5	12.60	78.3
		CTP	21.0 20.9	20.95	36.0
5	10.01.94	CTN	12.6 12.6	12.60	78.3
		CTP	20.9 20.8	20.85	36.0

10	Controls		Lot No.	Assigned Value
	CoagTrol N	CTN	537.001	73 - 99%
	CoagTrol P	CTP	541.034	29 - 39%

Table 34: Accelerated Stability

15

Assay: Prothrombin Time

Reagents: Innovin<sup>™</sup> PT Calibrator  
lot

20

Innovin<sup>™</sup> lot

Saline 0.9% lot

PILOT LOT 1TF - 12H1-86 H1 - 87

MLA E1000C Software version

5.00E P46

25

TEMPERATURE

INCUBATED: 20°C

30

Mean baseline clotting time (secs)	10.69	PT%	100
---------------------------------------	-------	-----	-----

Table 34 Continued

35

Clotting Time of neat plasma (seconds)						
Number of Day	Vial 1	Vial 2	Vial 3	Mean	% Change from Mean	PT %
5	11.2	11.3	11.3	11.10	3.84	92.8
	10.9	10.9	11.0			

51

Table 34 Continued

5	11	11.4	11.5	11.3	11.28	5.52	90.3
		11.3	11.1	11.1			
5	15	11.5	11.6	11.5	11.38	6.45	88.9
		11.3	11.2	11.2			
10	20	11.3	11.2	—	11.20	4.77	91.4
		11.2	11.1	—			
10	32	11.5	11.9	11.7	11.62	8.70	85.8
		11.5	11.5	11.6			
15	37	11.7	11.7	11.6	11.52	7.76	87.1
		11.4	11.4	11.3			
15	45	11.9	11.8	11.9	11.73	9.73	84.5
		11.7	11.5	11.6			
15	56	12.1	—	—	12.00	12.25	81.3
		11.9	—	—			

Failed stability after (days)	45
-------------------------------	----

## 20 Table 35: Accelerated Stability

Assay: Prothrombin Time

25 Reagents: Innovin™ PT Calibrator lot PILOT LOT 1  
 Innovin™ lot TFS - 12  
 Saline 0.9% lot H1-86 H1-87

MLA E1000C Software version 5.00E P46

30 TEMPERATURE INCUBATED: 30 °C

Mean baseline clotting time (secs)	10.69	PT%	100
------------------------------------	-------	-----	-----

35

52

Table 35 Continued

Clotting Time of neat plasma (seconds)							
Number of Day	Vial 1	Vial 2	Vial 3	Mean	% Change from Mean	PT %	
5	3	11.3	11.3	11.3	11.13	4.12	92.4
		10.9	11.0	11.0			
	4	11.7	11.6	11.5	11.47	7.30	87.7
		11.5	11.2	11.3			
	7	11.8	11.9	11.7	11.65	8.98	85.4
		11.5	11.5	11.5			
10	8	11.7	11.7	11.7	11.63	8.79	85.7
		11.7	11.5	11.5			
	9	11.9	11.9	11.8	11.73	9.73	84.5
		11.8	11.6	11.4			
	10	11.9	11.9	11.8	11.73	9.73	84.5
		11.6	11.8	11.4			
15	14	12.1	12.2	12.1	12.05	12.72	80.8
		12.0	12.0	11.9			
	15	12.6	12.8	—	12.65	18.33	74.6
		12.5	12.7	—			
	16	12.2	12.2	12.4	12.15	13.66	79.7
		12.0	12.1	12.0			
20	18	12.5	12.4	12.3	12.33	15.34	77.8
		12.3	12.3	12.2			
	20	12.6	12.8	—	12.60	17.87	75.1
		12.5	12.5	—			
25	Failed stability after (days) 10						

53

Table 36: Accelerated Stability

Assay: Prothrombin Time

5 Reagents: Innovin<sup>™</sup> PT Calibrator lot PILOT LOT 1  
 Innovin<sup>™</sup> lot TFS - 12  
 Saline 0.9% H1-86 H1 - 87

MLA E1000C Software version

5.00E P46

10

TEMPERATURE INCUBATED: 37°C

Mean baseline clotting time (secs)	10.69	PT%	100
------------------------------------	-------	-----	-----

15

Clotting Time of neat plasma (seconds)						
Number of Day	Vial 1	Vial 2	Vial 3	Mean	% Change from Mean	PT %
1	11.4	11.5	11.3	11.30	5.70	90.0
	11.2	11.2	11.2			
2	11.6	11.6	11.6	11.45	7.11	88.0
	11.3	11.3	11.3			
3	11.6	11.8	11.7	11.62	8.70	85.8
	11.6	11.5	11.5			
4	11.8	11.9	11.9	11.68	9.26	85.1
	11.7	11.7	11.6			
5	12.1	12.3	12.2	12.08	13.00	80.4
	12.0	12.0	11.9			
6	12.6	12.3	12.5	12.32	15.25	77.9
	12.2	12.0	12.3			

30

Failed stability after (days)	4
-------------------------------	---

54

Table 37: Accelerated Stability

Assay: Prothrombin Time

5 Reagents: Innovin<sup>™</sup> PT Calibrator lot PILOT LOT 1  
 Innovin<sup>™</sup> lot TFS - 12  
 Saline 0.9% lot H1-86 H1-87

10 MLA E1000C Software version 5.00E P46  
 TEMPERATURE INCUBATED: 50°C

Mean baseline clotting time (secs)	10.69	PT%	100
------------------------------------	-------	-----	-----

15

Clotting Time of neat plasma (seconds)						
Number of Hours	Vial 1	Vial 2	Vial 3	Mean	% Change From Mean	PT %
20 2	11.3	11.2	11.2	11.12	4.02	92.6
	11.1	11.0	10.9			
4	11.5	11.6	11.4	11.32	5.89	89.7
	11.2	11.1	11.1			
6	11.6	11.6	11.7	11.45	7.11	88.0
	11.2	11.2	11.4			
25 8	11.7	11.7	11.8	11.67	9.17	85.2
	11.6	11.5	11.7			

30

Failed stability after (hours)	8
--------------------------------	---



We claim:

1. A calibrator for a prothrombin time assay comprising:

- a) normal pool plasma selected from the group  
5 consisting of citrated plasma and citrate based anticoagulant plasma;
  - b) a quantity of a coagulation factor, which when added to said plasma, is sufficient to increase the %PT value of the plasma to about  
10 100% prior to lyophilization of said normal pool plasma and added coagulation factor.
2. The calibrator of claim 1 wherein the coagulation factor is selected from the group consisting of human rFVII, human rFVIIa, rFVII  
15 purified from at least one human plasma source, rFVIIa purified from at least one human plasma source, rFVII purified from any species' plasma source, rFVIIa purified from any species' plasma source, rFVIIa from any species' plasma source, rFVII from any species'  
20 plasma source, and any reagent with substantially the same functional activity of the aforementioned coagulation factor, including any dilution and mutation of any of such coagulation factors.
3. The calibrator of claim 2 wherein the  
25 increased %PT of said plasma and added coagulation factor is about 100% PT after lyophilization.

4. The calibrator of claim 3 for use with a PT reagent.

5. The calibrator of claim 3 for use with a recombinant tissue factor PT reagent.

5 6. The calibrator of claim 5 wherein the recombinant tissue factor PT reagent is selected from the group consisting of Innovin™ PT reagent and Ortho® RecomboPlastin™ PT reagent.

7. A method of preparing a calibrator for use in  
10 the prothrombin time assay comprising the steps of:

- a) collecting a normal pool of plasma selected from the group consisting of citrated plasma and citrate based anticoagulant plasma;
- b) adding a quantity of a coagulation factor to  
15 said plasma, which is sufficient to increase the %PT of said plasma and added coagulation factor to about 100% PT prior to lyophilization.

8. The method of claim 7 wherein the coagulation factor is selected from the group consisting of human  
20 rFVII, human rFVIIa, rFVII purified from at least one human plasma source, rFVIIa purified from at least one human plasma source, rFVII purified from any species' plasma source, rFVIIa purified from any species' plasma source, rFVIIa from any species' plasma source,  
25 rFVII from any species' plasma source, and any reagent with substantially the same functional activity of the

57

aforementioned coagulation factor, including any dilution and mutation of any of such coagulation factors.

9. The method of claim 8 wherein the increased  
5 %PT of said plasma and added coagulation factor is approximately 100% after lyophilization.

10. The method of claim 8 wherein the calibrator is used with Innovin™ PT reagent.

11. The method of claim 9 wherein the calibrator  
10 is used with a recombinant tissue factor PT reagent.

12. The method of claim 11 wherein the recombinant PT reagent is selected from the group consisting of Innovin™ PT reagent and Ortho® RecomboPlastin™ reagent.


15 13. A calibrator for a coagulation factor assay comprising:

- a. normal pool plasma selected from the group consisting of citrated plasma and citrate based anticoagulant plasma;
- 20 b. a quantity of coagulation factor, which when added to said plasma, is sufficient to increase the percentage of said coagulation factor of the plasma to about 100% prior to lyophilization of said normal pool plasma and added coagulation  
25 factor.

14. A method of preparing a calibrator for use  
in a coagulation factor assay comprising the steps of:

- 5       a. collecting a normal pool plasma selected from  
the group consisting of citrated plasma and  
citrate based anticoagulant plasma;
- 10       b. adding a quantity of coagulation factor,  
which when added to said plasma, is sufficient to  
increase the percentage of said coagulation  
factor of the plasma to about 100% prior to  
lyophilization of said normal pool plasma and  
added coagulation factor.

1/5

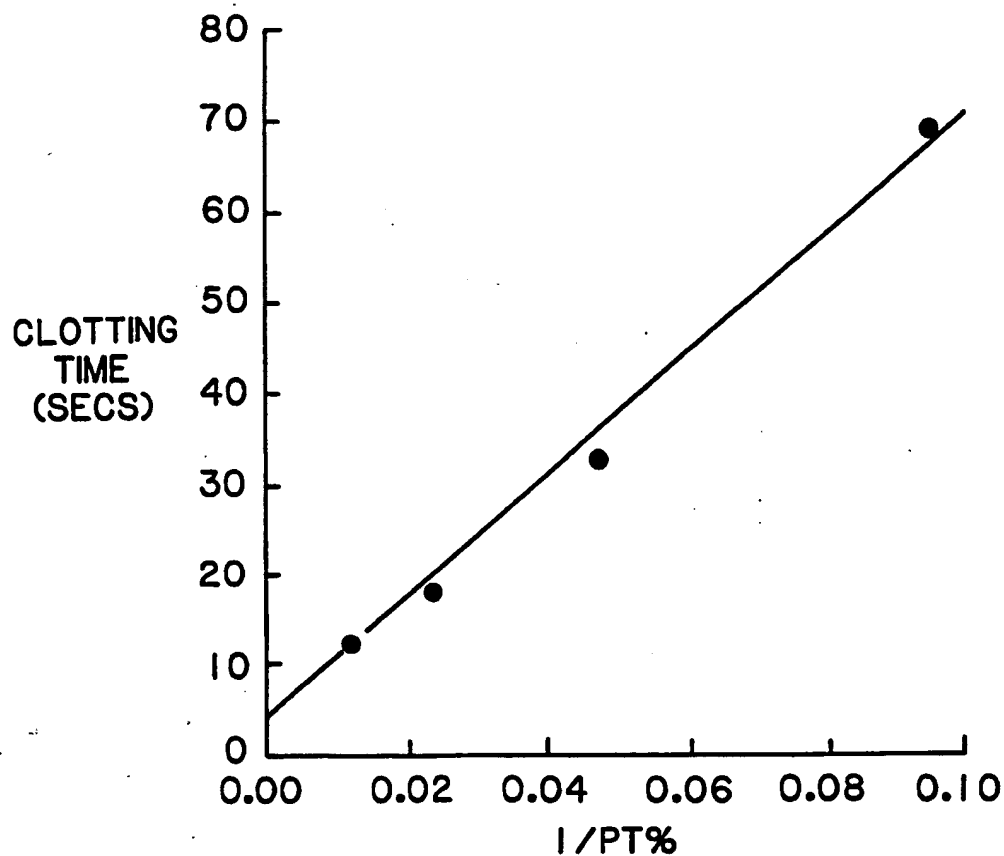


## DILUTION PT% CLOTING TIME

NEAT	85	12.2 secs
1 IN 2	42.5	18.0 secs
1 IN 4	21.25	32.7 secs
1 IN 8	10.625	69.0 secs

FIG. 1


2 / 5

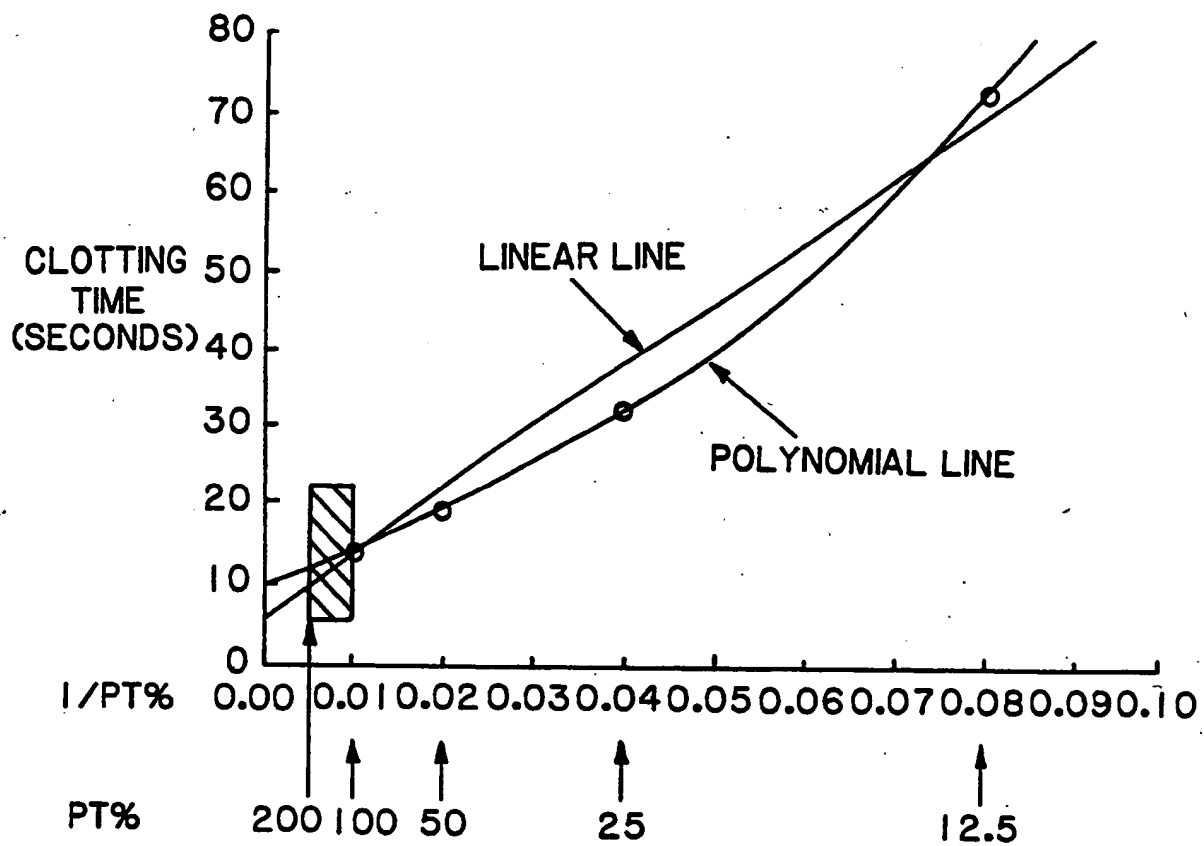

PT% VALUES FROM 100-200%

FIG. 2

3/5

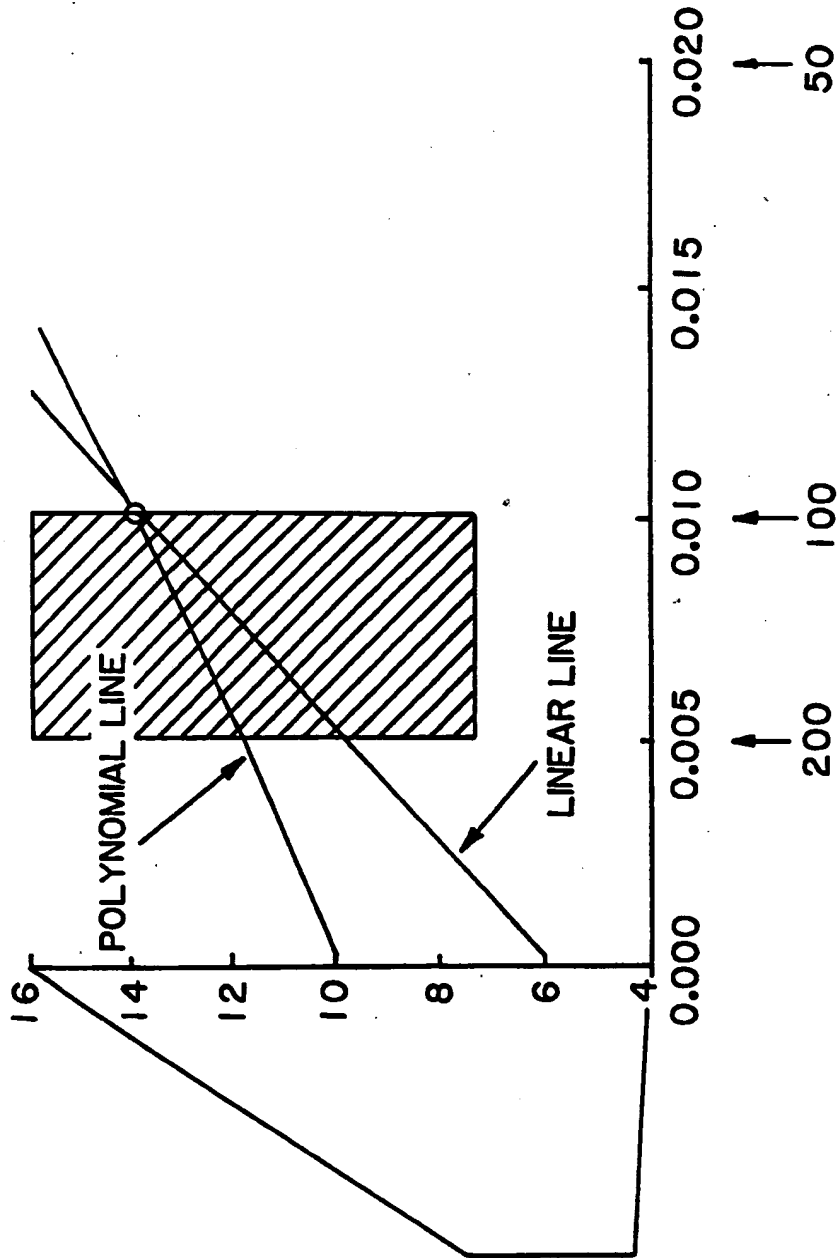

 PT% VALUES FROM 100-200%

FIG. 3

4 / 5

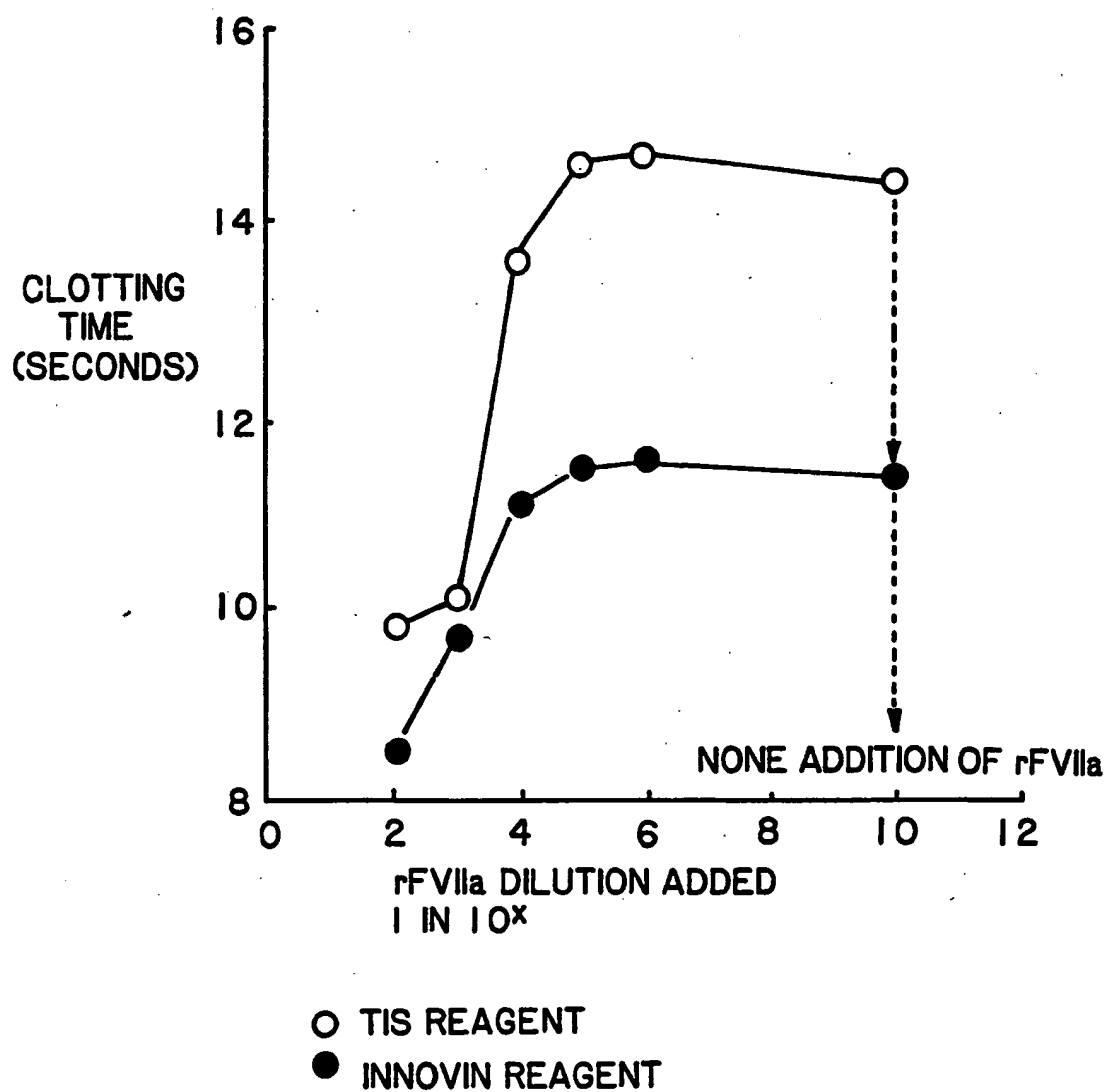


FIG. 4



5 / 5

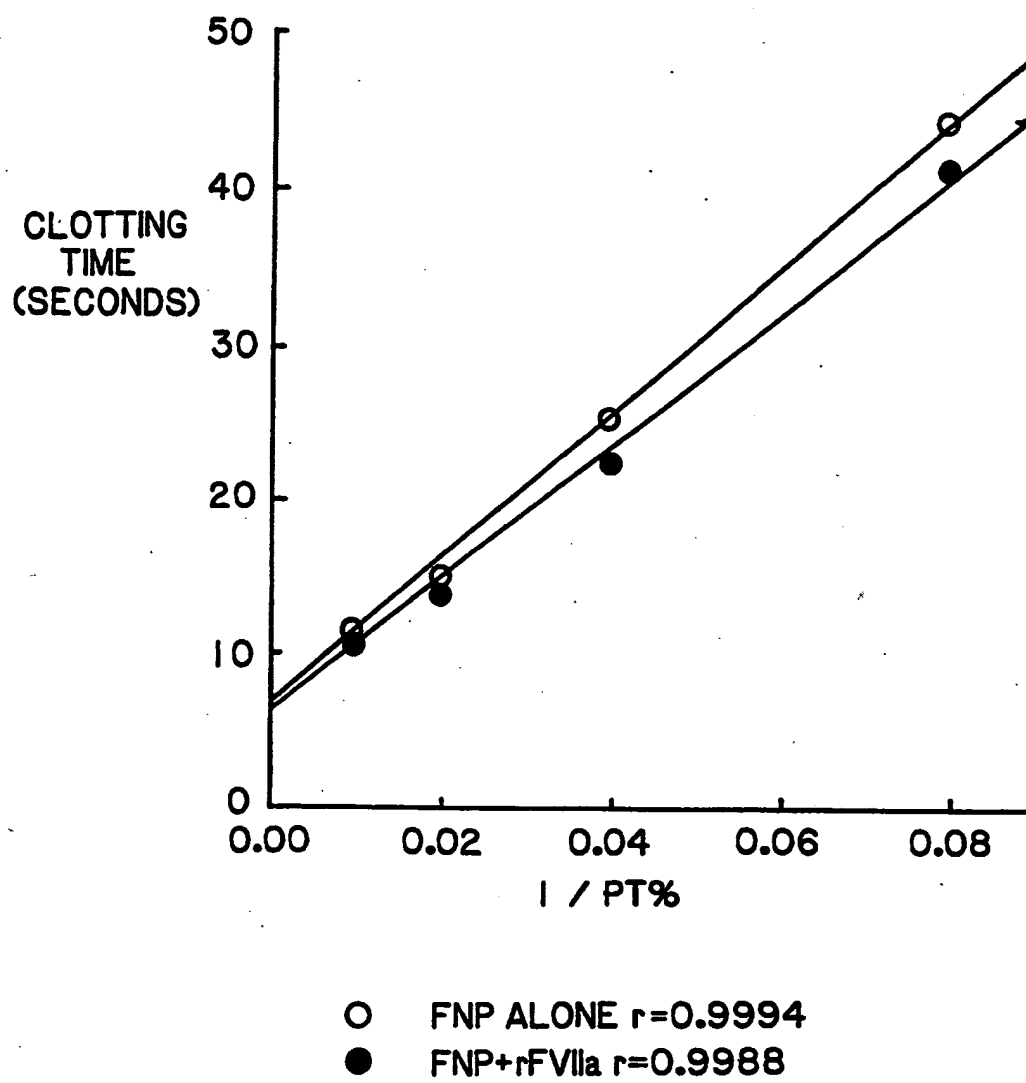


FIG. 5

# INTERNATIONAL SEARCH REPORT

Intern: J Application No

PCT/US 95/05195

**A. CLASSIFICATION OF SUBJECT MATTER**  
 IPC 6 G01N33/86 G01N33/96

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)

IPC 6 G01N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	EP-A-0 158 254 (BEHRINGWERKE AKTIENGESELLSCHAFT) 16 October 1985 ---	
P,X	KLIN. LAB., vol. 40, no. 7/8, August 1994 pages 619-628, U. SEYFERT ET AL. 'The determination of the Prothrombin Time in capillary blood using a Thromboplastin based on recombinant tissue factor and synthetic phospholipids.' see the whole document --- -/--	1-13

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

\* Special categories of cited documents :

- \*A\* document defining the general state of the art which is not considered to be of particular relevance
- \*E\* earlier document but published on or after the international filing date
- \*L\* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \*O\* document referring to an oral disclosure, use, exhibition or other means
- \*P\* document published prior to the international filing date but later than the priority date claimed

\*T\* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

\*X\* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

\*Y\* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

\*&\* document member of the same patent family

Date of the actual completion of the international search

23 August 1995

Date of mailing of the international search report

28.08.95

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
 NL - 2280 HV Rijswijk  
 Tel. (+ 31-70) 340-2040, Tx. 31 651 epo nl,  
 Fax: (+ 31-70) 340-3016

Authorized officer

Cartagena y Abella,P

# INTERNATIONAL SEARCH REPORT

Intern: 1 Application No

PCT/US 95/05195

## C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
P,A	<p>HÄMOSTASEOLOGIE,  vol. 14, no. 2, 5 July 1994  pages 90-99,  H. J. KOLDE ET AL. 'Erfahrungen mit einem  thromboplastin auf der basis von  rekombinantem gewebefaktor und  synthetischen phospholipiden.'  -----</p>	

**INTERNATIONAL SEARCH REPORT**

Information on patent family members

Int. Patent Application No

PCT/US 95/05195

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP-A-158254	16-10-85	DE-A- 3413311	17-10-85
		AU-B- 582570	06-04-89
		AU-A- 4089285	17-10-85
		CA-A- 1258221	08-08-89
		DE-A- 3563222	14-07-88
		JP-B- 7011524	08-02-95
		JP-A- 60230066	15-11-85
		US-A- 4784944	15-11-88
-----			

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**